

Gallium Artefacts on FIB-milled Silicon Samples


Joachim C. Reiner, Philipp M. Nellen, Urs Sennhauser

EMPA - Electronics/Metrology Laboratory, Ueberlandstr. 129, CH-8600 Dübendorf, Switzerland

Purpose


The Focused Ion Beam Machine (FIB) allows sample milling on the sub-micrometer scale by local sputtering. Here, for the case of silicon semiconductor sample preparation, two kinds of re-deposition artefacts are analysed, which can be disturbing particularly in failure analysis applications.

Analysis of Artefacts


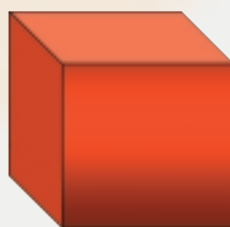


amount of Ga for 1 lamella:
 30 min, 20 nA
 $\sim 2.55 \cdot 10^{-8}$ g Ga
 $\sim 4400 \mu\text{m}^3$
 $\sim (16.4 \mu\text{m})^3$

TEM lamella preparation by FIB

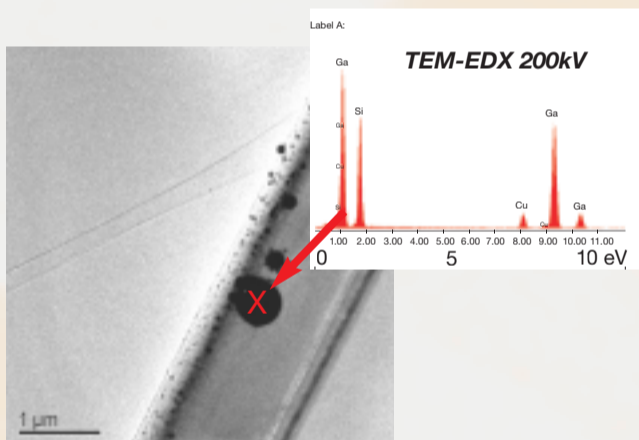

TEM bright-field image with preparation artefacts
 EDX spectrum of large artefact

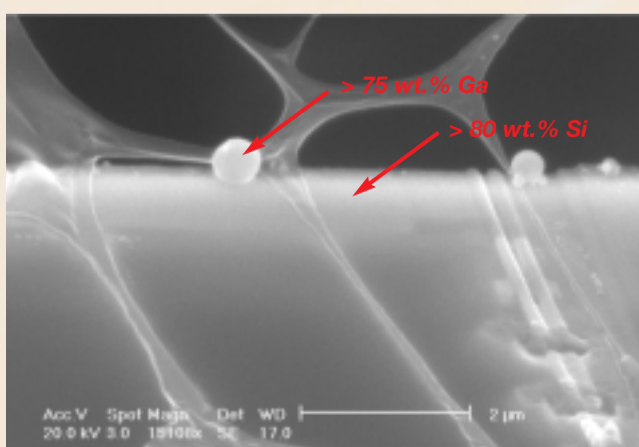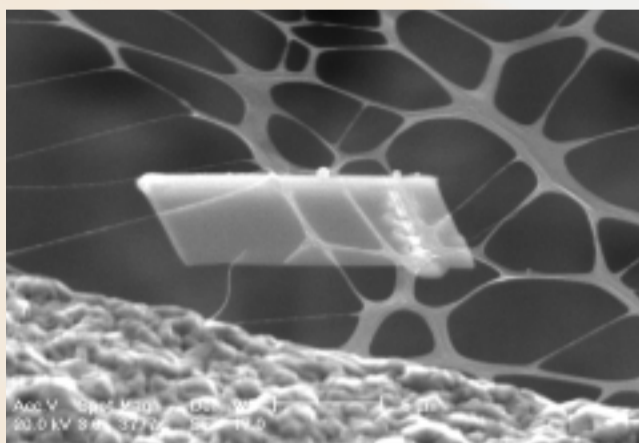

Sample lamella on TEM carrier (SEM, EDX) with Ga droplets

Experiments

Experiment 1: Ga-Beam-Induced Re-Deposition

Results:

1. A closed film re-deposition appears, expected to consist mainly of silicon.
2. Ga droplets form.

Experiment 2: Prevention of Re-Deposition

Result:

3. No re-deposition is observed in the shade.

Summary

- Two kinds of re-deposition: film & droplets.
- Preparation artefact «Droplets» consist of > 75 wt.% Ga. Such Ga droplets may mislead in failure analysis.
- Re-deposition is induced by the Ga-beam.
- Re-deposition prevention by keeping point of interest in the shade.

Acknowledgement

We gratefully acknowledge the FIB, TEM and SEM work by Philippe Gasser, Giovanni Nicoletti and Stephan Meier, and helpful discussions with Peter Jacob, EMPA, and Prof. John Melngailis, University of Maryland, USA.